

1000W-2000W

4 Wavelength

808nm Diode Laser Machine


- ✓ Double Touch Screen Controller
- ✓ Painted ABS Housings
- ✓ 15.6 Inch Color Touch Screen

755nm+808nm+940nm+1064nm



12*28mm



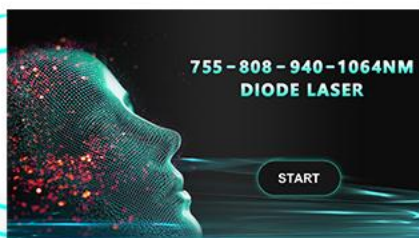

Large Screen

15.6 Inch Color Touch Screen



Convenient
Operation



Fastresponse
Time





Intelligent
Processing System

Color LCD Display

Color LCD Display Light On The Handle
Status And Parameters
Touchable Display





4 Wavelengths

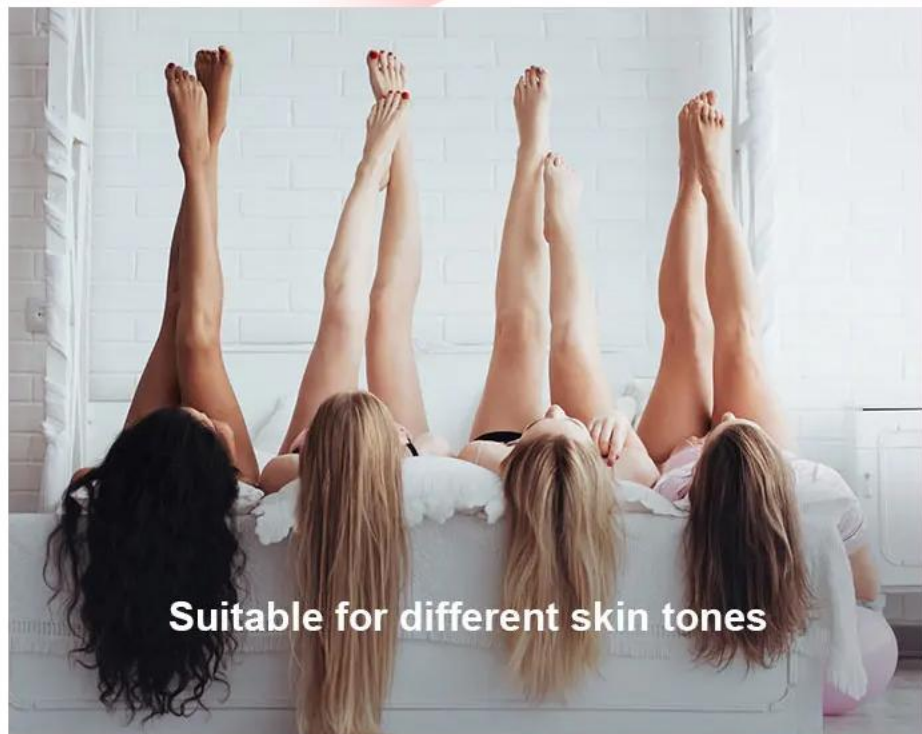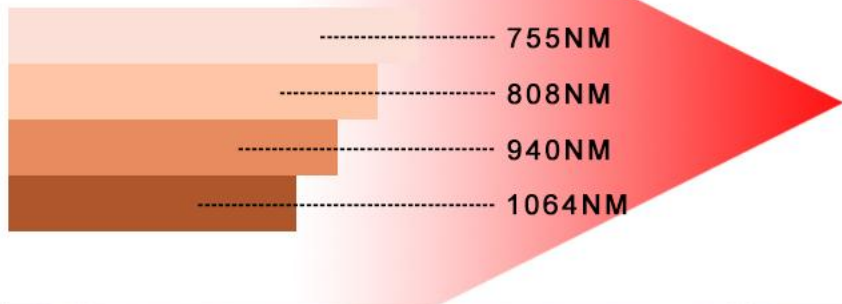

4 Wavelength Diagram

The Gold Standard For Hair Removal Treatment
Is Safe For The Epidermis

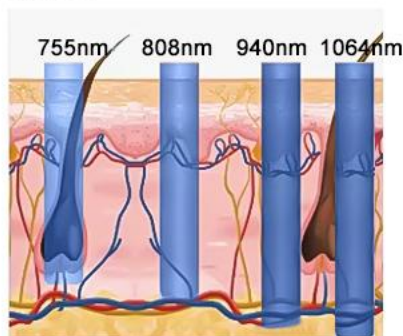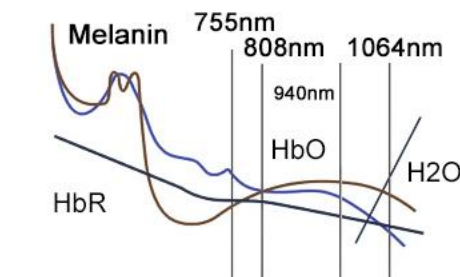

Why Choose Us

Refrigeration System



Double TEC

Imported From Japan, It Is Safer And Can control The Water Temperature At 33 Degrees within 5/6 Hours To Ensure That The Machine Can Run Continuously Within 24 Hours.

TEC+Water+Air cooling+Semiconductor refrigeration, All-Round Guarantee That The machine Always Starts At The Right Temperature, Ensuring A Comfortable Treatment Experience.

10-15 Years Longer Life.

Compressor

Imported From Japan, Used For Machines above 2000W, The Cooling Effect Is Better 1.5 Times That Of TEC.

500W Dc Double-Bearing Compressor With a Service Life Of 5 Years.

Small Size, High Energy, Low Fluorine Content. Shipping Is Not Affected.



Compared With Ordinary AC Compressors, It Is More Environmentally.



100 Million Shots Warranty




Light Beam Focusing Technology

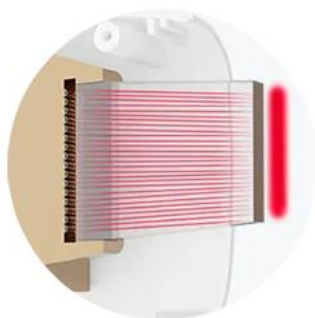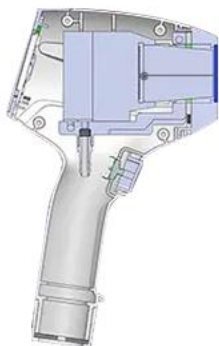

Laser Bar + German
Asphericoptical Microlens.

DPL

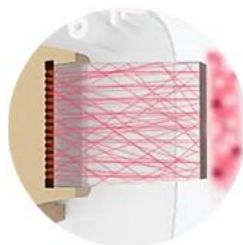


Just The Laser Bar.

OTHER



Aspheric microlens beam collimation technology the beam is collimated by aspherical microlenses, providing higher beam quality.

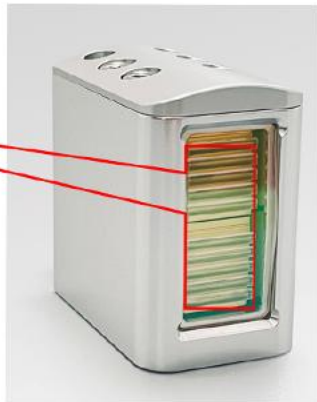

New Italy imported FLuid-o-tech water pump comply with international standards, strictly tested & quality control, with better cooling system, quiet working, strong motivation and excellent durability.

Uses Ion double filter to filter metal ions to ensure good water circulation, avoid internal laser corrosion and extend the service life of the instrument.



The LED Screen Displays

Healingenergy And The Number Of Rounds Used, For Easy Viewing At Any Time



Sapphire Crystal

Best Light And Cold Conductivity, Ensuring The Effectiveness And Comfort Of The Treatment

Cooling Effect

The Handle Is Frosted For One Minute.
Making Your Treatment More Comfortable And Painless

Just One Minute!

Quick Cooling Down To

-6°C



For White Skin

(Fine, Golden Hair Removal)



For Yellow Neutral Skin

(Brown Hair Removal)



For Black Skin

(Black Hair Removal)

Common Body Areas Covered Under Diode Laser Hair Removal



Treatment Effect



Working Interface

- 1、Start Interface; After Booting, The Startup Interface Is Displayed, As Shown Below..
- 2、Select The Corresponding Gender, Skin Color, Region, And Hair Color, And Click The "Enter" Key To Enter The Main Operation Interface.

